



VMRD

PO Box 502, Pullman, WA 99163 USA

Telephone: 509-334-5815

Fax: 509-332-5356

E-mail: vmrd@vmrd.com

Web site: <http://www.vmrd.com>

Certificate of Analysis

FELINE INFECTIOUS PERITONITIS (FIP)

Direct FA Conjugate


CATALOG NO.: 210-31-FIP

VOLUME: 10 ml

LOT: P081031-001

EXPIRATION: 3 November 2010

AGENT: Feline Infectious Peritonitis (FIP)

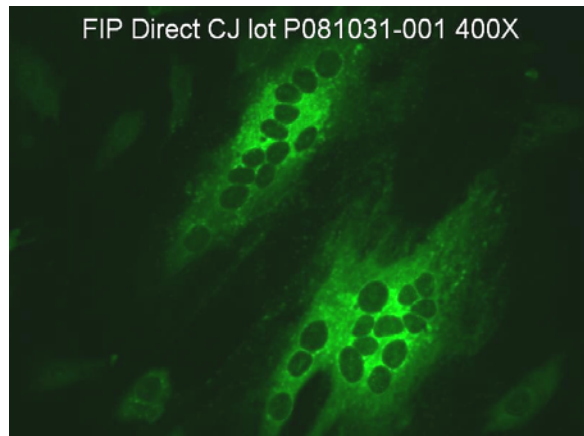

DESCRIPTION: Anti-FIP polyclonal antiserum conjugated to fluorescein isothiocyanate. Porcine and feline origin. Multivalent, staining both types 1 and 2. Ready to use.

QUALITY CONTROL METHOD: Direct FA on infected cell cultures using VMRD, Inc. FIP-1 (catalog no. SLD-FAC-FIP1) and FIP-2 (catalog no. SLD-FAC-FIP2) control slides to detect binding.

Specific Reaction: 3-4+ on FIP1 and FIP2 positive cells; negative on negative cells.

Other Comment: No background.

PATTERN OF FLUORESCENCE: Diffuse granular cytoplasmic.

INTENDED USE: Suitable for detecting FIP in tissues or cell cultures. Cross-reacts with Transmissible Gastroenteritis (TGE) and Canine Coronavirus (CCV).

STORAGE: This conjugate is provided in liquid form and should be stored at 2-7°C. **DO NOT FREEZE!** If conjugate becomes cloudy, it should be discarded. This conjugate contains 0.09% sodium azide as a preservative.

FOR IN VITRO LABORATORY USE ONLY.

WARRANTY: VMRD, Inc. warrants that this product is as described in the quantity and contents stated on the label at the time of delivery to the customer. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE MADE BEYOND THE LABEL DESCRIPTION, INCLUDING WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR USE. Remedy is limited to replacement of the product or refund of the purchase price. VMRD, Inc. is not liable for property damage, personal injury, or economic loss caused by the product. The information listed in this information sheet is provided for reference only, and should not be substituted for the user's own incoming material quality control.

RECOMMENDED STAINING PROCEDURE FOR DIRECT FA:

1. Air dry smears or tissue sections for at least 30 minutes at room temperature (do not dry cell cultures!).
2. Fix smears or tissue sections on slides for 20 minutes in acetone-methanol (75/25) at room temperature. Cell cultures should be rinsed with PBS and fixed in pure acetone for 10 minutes at room temperature. After fixation and before staining, slides should be dried for 10 minutes in a dry 37°C incubator.
3. Stain slides with 50-75 µl conjugate for 30 minutes at 37°C in humid chamber.
4. Gently rinse slides briefly in FA Rinse Buffer, pH 9.0 (VMRD catalog no. 210-90-RB) and then soak for 10 minutes in FA Rinse Buffer, pH 9.0.
5. Drain slides and dry back and edges with a paper towel. Do not allow stained surface to dry. Do not rinse with water.
6. Mount with FA Mounting Fluid [glycerol/FA rinse buffer, pH 9.0, (50/50)] (VMRD catalog no. 210-92-MF) and view with good quality fluorescence microscope at 100X-250X. Confirmation may be made at 400X.

PHOSPHATE BUFFERED SALINE (PBS) SOLUTION (pH 7.2):

- Na₂HPO₄1.19 gm
- NaH₂PO₄0.22 gm
- NaCl8.55 gm
- DI/dH₂OQ.S. to 1 liter

4X FA RINSE BUFFER (pH 9.0):

- Na₂CO₃11.4 gm
- NaHCO₃33.6 gm
- NaCl8.5 gm
- DI/dH₂OQ.S. to 1 liter

Final pH should be 9.0-9.5. This is a 4X concentrate and should be diluted 1:4 with DI/distilled water for use as a working buffer. Keep in a tightly stoppered container at room temperature. FA MOUNTING FLUID is made by mixing glycerol and FA Rinse Buffer, pH 9.0, in equal proportions.